

AseptiQuik G®

Series Connectors

AseptiQuik® G Connectors enable quick and easy intrinsic sterile connections, even in non-sterile environments. The easy-to-use genderless design simplifies system integration and minimizes the risk of operator error. The connectors' robust construction provides enhanced user confidence and reliable performance without the need for clamps, fixtures or tube welders. Biopharmaceutical manufacturers can maintain their closed systems using intrinsic sterile connectors with the quality and market availability expected from the leader in single-use connections.

SPECIFICATIONS

OPERATING PRESSURE

Up to 60 psi, 4.1 bar
Up to 75 psi, 5.1 bar for 48 hours (excluding AQGMPC)

OPERATING TEMPERATURE

34°F to 104°F (1°C to 40°C)

TERMINATIONS

1/4", 3/8", 1/2", 3/4" ID hose barb (6.4 mm, 9.5 mm, 12.7 mm, 19.0 mm) and 3/4", 1-1/2" sanitary

MATERIALS

Main Components: White polycarbonate (PC) for standard and HT versions
Off-white polyphenylsulfone (PPSU) for standard and HT versions; polycarbonate weld-plate (non-flow path material) for PPSU HT

Pull Tabs/Caps: Polycarbonate:

Blue = standard temp, White = high temp,
Purple = standard temp PPSU, Gold = high temp PPSU

Seals: Silicone (clear), platinum-cured

Membranes: Polyethylene (standard PC and PPSU versions) Hydrophobic polyethersulfone (HT PC and PPSU versions) PTFE slider sticker (HT PC and PPSU versions)

STERILIZATION

Standard (blue) and PPSU (purple) Versions:


Gamma: up to 50kGy

High Temperature (white) and PPSU (gold) Versions:

Gamma: up to 50kGy

Autoclave: One cycle up to 266°F (130°C) for 60 minutes

WARNING: Pressure, temperature, chemicals, and operating environment can affect the performance of couplings. It is the customer's responsibility to test the suitability of CPC's products in their own application conditions.

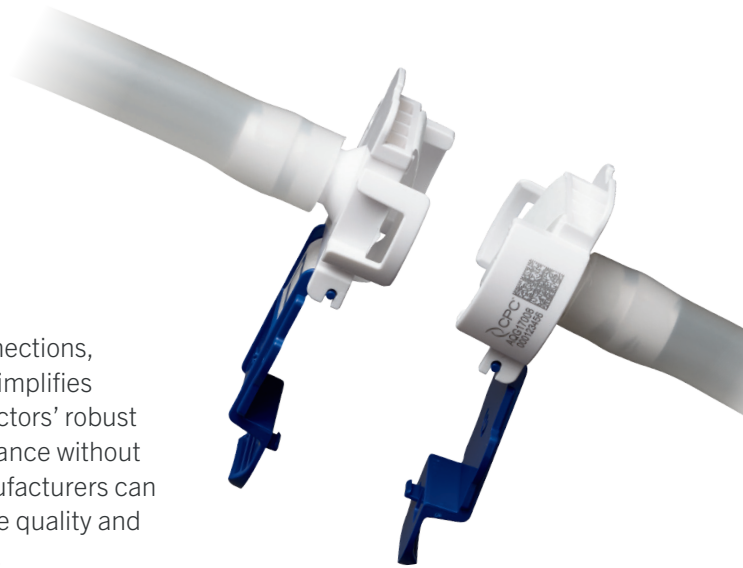

FEATURES	BENEFITS
✓ Genderless	Eases single-use systems specifications with one-part number for both halves
✓ FLIP-CLICK-PULL	Intuitive three-step connection process reduces risk of operator error
✓ CPC Click	Audible confirmation of connection with no additional hardware required
✓ Chemical Compatibility and pH Range	AseptiQuik PPSU enables genderless connection for a greater range of chemical applications, offering versatile connections across downstream processes with a pH range from 2 to 12
✓ BPA-free	AseptiQuik PPSU meets a broader range of Application Requirements


TYPICAL FLOW RATE

Cv Value Range: 1.5 - 31

Cv values represent the approximate expected flow rate in gallons per minute of water at room temperature for a 1 PSI pressure drop. The flow is generally constrained by the smallest diameter, which in some cases will be the termination diameter and not the Nominal Flow Path.

ASEPTIQUIK G® SERIES DIMENSIONS

PRODUCT DIMENSIONS

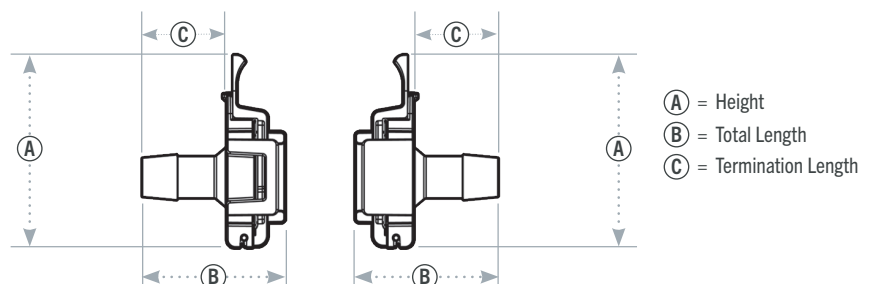

POLYCARBONATE with blue pull tabs For gamma irradiation applications.

AseptiQuik G® Series Connectors



TERMINATION	METRIC EQ.	PART NO.	Ⓐ HEIGHT	Ⓑ TOTAL LENGTH	Ⓒ TERMINATION LENGTH
1/4" HOSE BARB	6.4 mm	AQG17004	2.62" (66.6 mm)	1.62" (41.2 mm)	0.82" (20.8 mm)
3/8" HOSE BARB	9.5 mm	AQG17006	2.62" (66.6 mm)	1.62" (41.2 mm)	0.82" (20.8 mm)
1/2" HOSE BARB	12.7 mm	AQG17008	2.62" (66.6 mm)	1.95" (49.5 mm)	1.15" (29.2 mm)
3/4" HOSE BARB	19.1 mm	AQG17012	2.62" (66.6 mm)	2.36" (59.9 mm)	1.56" (39.6 mm)
3/4" SANITARY	19.1 mm	AQG33012	2.62" (66.6 mm)	1.66" (42.2 mm)	0.86" (21.8 mm)
1-1/2" SANITARY	38.1 mm	AQG33024	2.62" (66.6 mm)	1.93" (48.9 mm)	1.13" (28.6 mm)
MPC INSERT		AQG17MPC	2.62" (66.6 mm)	1.50" (38.1 mm)	0.70" (17.8 mm)

POLYCARBONATE with white pull tabs For autoclave or gamma irradiation applications.



TERMINATION	METRIC EQ.	PART NO.	Ⓐ HEIGHT	Ⓑ TOTAL LENGTH	Ⓒ TERMINATION LENGTH
1/4" HOSE BARB	6.4 mm	AQG17004HT	2.62" (66.6 mm)	1.62" (41.2 mm)	0.82" (20.8 mm)
3/8" HOSE BARB	9.5 mm	AQG17006HT	2.62" (66.6 mm)	1.62" (41.2 mm)	0.82" (20.8 mm)
1/2" HOSE BARB	12.7 mm	AQG17008HT	2.62" (66.6 mm)	1.95" (49.5 mm)	1.15" (29.2 mm)
3/4" HOSE BARB	19.1 mm	AQG17012HT	2.62" (66.6 mm)	2.36" (59.9 mm)	1.56" (39.6 mm)
3/4" SANITARY	19.1 mm	AQG33012HT	2.62" (66.6 mm)	1.66" (42.2 mm)	0.86" (21.8 mm)
1-1/2" SANITARY	38.1 mm	AQG33024HT	2.62" (66.6 mm)	1.93" (48.9 mm)	1.13" (28.6 mm)
MPC INSERT		AQG17MPCHT	2.62" (66.6 mm)	1.50" (38.1 mm)	0.70" (17.8 mm)

POLYPHENYLSULFONE with purple pull tabs For gamma irradiation applications.



TERMINATION	METRIC EQ.	PART NO.	Ⓐ HEIGHT	Ⓑ TOTAL LENGTH	Ⓒ TERMINATION LENGTH
1/4" HOSE BARB	6.4 mm	AQG17104	2.62" (66.6 mm)	1.62" (41.2 mm)	0.82" (20.8 mm)
3/8" HOSE BARB	9.5 mm	AQG17106	2.62" (66.6 mm)	1.62" (41.2 mm)	0.82" (20.8 mm)
1/2" HOSE BARB	12.7 mm	AQG17108	2.62" (66.6 mm)	1.95" (49.5 mm)	1.15" (29.2 mm)
3/4" HOSE BARB	19.1 mm	AQG17112	2.62" (66.6 mm)	2.36" (59.9 mm)	1.56" (39.6 mm)
3/4" SANITARY	19.1 mm	AQG33112	2.62" (66.6 mm)	1.66" (42.2 mm)	0.86" (21.8 mm)
1-1/2" SANITARY	38.1 mm	AQG33124	2.62" (66.6 mm)	1.93" (48.9 mm)	1.13" (28.6 mm)

POLYPHENYLSULFONE with gold pull tabs For autoclave or gamma irradiation applications.



TERMINATION	METRIC EQ.	PART NO.	Ⓐ HEIGHT	Ⓑ TOTAL LENGTH	Ⓒ TERMINATION LENGTH
1/4" HOSE BARB	6.4 mm	AQG17104HT	2.62" (66.6 mm)	1.62" (41.2 mm)	0.82" (20.8 mm)
3/8" HOSE BARB	9.5 mm	AQG17106HT	2.62" (66.6 mm)	1.62" (41.2 mm)	0.82" (20.8 mm)
1/2" HOSE BARB	12.7 mm	AQG17108HT	2.62" (66.6 mm)	1.95" (49.5 mm)	1.15" (29.2 mm)
3/4" HOSE BARB	19.1 mm	AQG17112HT	2.62" (66.6 mm)	2.36" (59.9 mm)	1.56" (39.6 mm)
3/4" SANITARY	19.1 mm	AQG33112HT	2.62" (66.6 mm)	1.66" (42.2 mm)	0.86" (21.8 mm)
1-1/2" SANITARY	38.1 mm	AQG33124HT	2.62" (66.6 mm)	1.93" (48.9 mm)	1.13" (28.6 mm)



[View Product Page](#)

DESIGN FOR LIFE



WARNING: The customer is solely responsible for determining the necessary or most suitable product or service for any particular use case or application. All such risks are assumed by the customer. Specifications subject to change without notice.
WARRANTY: For product warranty information, please contact: info@cpc-bio.com or +1 800-519-7633.
 For trademark information, please visit: www.cpc-bio.com/trademarks-and-patents
 © 2026 CPC BIOTECH. ALL RIGHTS RESERVED.
 REV 6/26

CPC BIOTECH
 2820 Cleveland Ave. N.
 Roseville, MN 55113
 Phone: +1 800-519-7633
 Email: info@cpc-bio.com

cpc-bio.com